# beurer


## SC 30

### Green Planet

|           |                                                                 |    |
|-----------|-----------------------------------------------------------------|----|
| <b>DE</b> | Elektrische Zahnbürste<br>Gebrauchsanweisung .....              | 6  |
| <b>EN</b> | Electric toothbrush<br>Instructions for use.....                | 13 |
| <b>FR</b> | Brosse à dents électrique<br>Mode d'emploi.....                 | 19 |
| <b>ES</b> | Cepillo de dientes eléctrico<br>Manual de instrucciones .....   | 25 |
| <b>IT</b> | Spazzolino elettrico<br>Istruzioni per l'uso .....              | 32 |
| <b>TR</b> | Elektrikli diş fırçası<br>Kullanım kılavuzu .....               | 38 |
| <b>RU</b> | Электрическая зубная щетка<br>Инструкция по<br>применению ..... | 44 |
| <b>PL</b> | Elektryczna szczoteczka<br>do zębów<br>Instrukcja obsługi ..... | 51 |
| <b>NL</b> | Elektrische tandenborstel<br>Gebruiksaanwijzing .....           | 58 |
| <b>DA</b> | El-tandbørste<br>Betjeningsvejledning.....                      | 64 |
| <b>SV</b> | Eltandborste<br>Bruksanvisning.....                             | 70 |
| <b>NO</b> | Elektrisk tannbørste<br>Bruksanvisning.....                     | 76 |
| <b>FI</b> | Sähköhammasharja<br>Käyttöohje.....                             | 82 |

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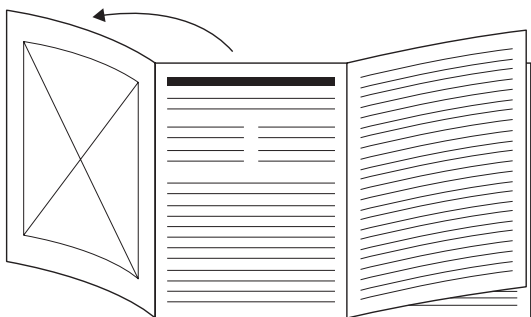Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the next user with these instructions for use.

## CONTENTS

|                                           |    |
|-------------------------------------------|----|
| 1. Warnings and safety notes .....        | 13 |
| 2. Included in delivery .....             | 14 |
| 3. Signs and symbols .....                | 15 |
| 4. Intended purpose .....                 | 16 |
| 5. Device description.....                | 16 |
| 6. Usage .....                            | 16 |
| 6.1 Preparing for use.....                | 16 |
| 6.2 Starting use.....                     | 17 |
| 6.3 Travel lock .....                     | 17 |
| 7. Cleaning and maintenance.....          | 17 |
| 8. Accessories and replacement parts..... | 17 |
| 9. Disposal.....                          | 18 |
| 10. Technical specifications .....        | 18 |
| 11. Guarantee.....                        | 18 |

## 1. WARNINGS AND SAFETY NOTES

### **⚠ WARNING**

- Keep the device and the charging cable away from water or other liquids.
- Do not place the device near a sink or bathtub full of water.
- This device may be used by children over the age of 8 years and by people with reduced physical, sensory or mental capabilities or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the risks of using the device.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.

- Do not use the device if it shows signs of damage or does not function properly. In these cases, contact Customer Services.
- If the mains connection cable is broken, you must dispose of the charging cable.
- Only use the USB-C charging cable supplied with your device.
- Never use the device with the mains adapter near a bath, sink, shower or other containers filled with water or other liquids – risk of electric shock.
- The device may only be charged with a mains adapter (not included in delivery) with a safety extra-low voltage output of 5 V and max. 1 A.
- Under no circumstances should you open or repair the device yourself, otherwise proper functionality can no longer be guaranteed. Failure to comply will invalidate the warranty.

### Notes on handling batteries

#### **⚠ WARNING**

- If fluid from a battery comes into contact with your skin or eyes, wash the affected areas with water and seek medical assistance.
- Risk of explosion! Do not throw batteries into a fire.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Do not disassemble, split or crush the batteries.

#### **⚠ CAUTION**

- Protect batteries from excessive heat.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use (see Chapter “Preparing for use”).
- In order to achieve as long a battery service life as possible, fully charge the battery at least 2 times a year.

## 2. INCLUDED IN DELIVERY

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been

removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 electric toothbrush
- 2 brush heads
- 1 USB-C charging cable
- Instructions for use

### 3. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

|                                                                                                                          |                                                                                                                                       |
|--------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------|
|                                          | <b>DANGER</b>                                                                                                                         |
| Indicates an imminent danger. If it is not avoided, death or serious injury will occur.                                  |                                                                                                                                       |
|                                          | <b>WARNING</b>                                                                                                                        |
| Indicates a potentially imminent hazard. If it is not avoided, there is a risk of death or serious injury.               |                                                                                                                                       |
|                                          | <b>CAUTION</b>                                                                                                                        |
| Indicates a potentially imminent hazard. If it is not avoided, slight or minor injuries may result.                      |                                                                                                                                       |
|                                          | <b>NOTICE</b>                                                                                                                         |
| Indicates a potentially harmful situation. If it is not avoided, the device or something in its vicinity may be damaged. |                                                                                                                                       |
|                                        | <b>Product information</b><br>Note on important information                                                                           |
|                                        | Read the instructions                                                                                                                 |
|                                        | Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).                             |
|                                        | Do not dispose of batteries containing hazardous substances with household waste                                                      |
|                                        | The products demonstrably meet the requirements of the Technical Regulations of the EAEU.                                             |
|                                        | Standard for independent certification of recycled material                                                                           |
|                                        | <b>CE labelling</b><br>This product satisfies the requirements of the applicable European and national directives.                    |
|                                        | Separate the product and packaging elements and dispose of them in accordance with local regulations.                                 |
|                                        | Marking to identify the packaging material.<br>A = Material code, B = Material number:<br>1–7 = Plastics, 20–22 = Paper and cardboard |
|                                        | Conformity assessed marking for the United Kingdom                                                                                    |
|                                        | Protection class III device                                                                                                           |

|                                                                                  |                                                                                                             |
|----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------|
|   | Importer symbol                                                                                             |
|  | Manufacturer                                                                                                |
|  | Protection against temporary submersion. Device can be immersed to a depth of 1 metre for up to 30 minutes. |
|  | USB / Micro-USB connection                                                                                  |
|  | For indoor use only                                                                                         |
|  | Direct current                                                                                              |
|  | Packaging made from recycled material                                                                       |

## 4. INTENDED PURPOSE

Only use the device to clean teeth on humans. Do not use the device on animals! The device is not suitable for children under 8 years. The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

## 5. DEVICE DESCRIPTION

The associated drawings are shown on page 3.

- 1** Brush head
- 2** ON/OFF/Mode button
- 3** LED display
  - 4 modes: Clean, White, Polish, Gum Care
  - Battery level
- 4** USB-C charging cable

## 6. USAGE

### 6.1 Preparing for use

#### Charging the electric toothbrush **A**

Fully charge the electric toothbrush before first use. To do so, proceed as follows:

1. Connect the USB-C charging cable to the electric toothbrush. Plug the USB-C charging cable into a USB port or into a suitable mains adapter and plug it into a socket.
2. Whilst charging, the LED pulsates red. As soon as the electric toothbrush is fully charged, the LED lights up green continuously.

#### Fitting the brush head **B**

Place the brush head **1** onto the device from above until it clicks into place. Make sure that the brush head **1** points in the same direction as the ON/OFF/Mode button **2**.

#### Indicator bristles

Indicator bristles lose their colour over time depending on how each person uses the device, thus indicating the best time to change the brush head **1**.

We recommend changing the brush head **1** after 3 months. If the bristles of the brush head are worn or frayed, they may need to be replaced more frequently.

## 6.2 Starting use

1. To switch on the electric toothbrush, briefly press the ON/OFF/Mode button **2**. The electric toothbrush pauses briefly every 30 seconds to remind you to clean all 4 jaw quadrants for the same amount of time **C**. The electric toothbrush has an automatic timer and vibrates two times after 2 minutes. The toothbrush switches off automatically after 6 minutes.
2. Select your desired mode by briefly pressing the ON/OFF/Mode button **2**. The currently selected mode is shown accordingly by the corresponding mode LED. The following modes are available:

|          |                                 |
|----------|---------------------------------|
| Clean    | For daily cleaning              |
| White    | For whiter teeth                |
| Polish   | For removing superficial stains |
| Gum Care | For looking after the gums      |

3. To cancel the application prematurely, briefly press the ON/OFF/Mode button **2**.

## 6.3 Travel lock

The travel lock prevents the device from being switched on by accident (e.g. on trips, when being carried in a suitcase).

1. To activate the travel lock, press and hold the ON/OFF/Mode button **2** for 5 seconds. During activation, the LED flashes green 5 times.
2. If you press the ON/OFF/Mode button **2** briefly when the device is locked, the LED will also flash green.
3. To deactivate the travel lock, press and hold the ON/OFF/Mode button **2** again for 5 seconds. During deactivation, the LED lights up green continuously. The LED goes out as soon as the travel lock has been deactivated. The device starts up.

## 7. CLEANING AND MAINTENANCE

Clean the device and the brush head **1** after each use. To do so, proceed as follows:

1. Pull the brush head **1** upwards and off.
2. Clean the brush head **1** under running water. Make sure to clean the inside of the brush head **1** as well. **D**
3. Clean the device under running water. However, do not submerge the device in water. Water might enter the device and damage it. **E**
4. Regularly clean the charger **5** with a damp cloth and a mild soap solution. Do not use any abrasive cleaning products. Keep the inductive charger **5** and the USB charging cable with mains adapter away from water.

## 8. ACCESSORIES AND REPLACEMENT PARTS

To purchase accessories and replacement parts, please visit [www.beurer.com](http://www.beurer.com) or contact the corresponding service address (according to the service address list) in your country. Accessories and replacement parts are also available from retailers.

| Designation                                | Item number and/or order number |
|--------------------------------------------|---------------------------------|
| Classic Clean brush heads 4-pack           | 105.21                          |
| Premium Whiten brush heads 4-pack          | 105.31                          |
| Premium Plaque Eraser brush heads 4-pack   | 105.22                          |
| Premium Protection Plus brush heads 4-pack | 105.23                          |

9. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Observe the local regulations for material disposal. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal. You can obtain the location of collection points for old devices from the local authorities.

The empty, completely flat batteries must be disposed of in specially designated collection boxes, at recycling points or at electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on batteries containing harmful substances:

Pb = battery contains lead

Cd = battery contains cadmium

Hg = battery contains mercury

Removing the battery

- 1. Use a screwdriver to open the electric toothbrush. Insert the screwdriver into the USB-C port.
- 2. Now apply a little pressure and turn anti-clockwise. The end cap on the toothbrush can now be removed.
- 3. Fold the bottom of the toothbrush to the side.
- 4. Pull out the inner housing of the toothbrush.
- 5. Remove the battery from the inner housing.
- 6. Dispose of the battery separately in specially designated collection boxes, at recycling points or at electronics retailers.

10. TECHNICAL SPECIFICATIONS

|                                                             |                                            |
|-------------------------------------------------------------|--------------------------------------------|
| Type                                                        | SC 30                                      |
| Dimensions                                                  | Device: 183 x 29 mm                        |
| Weight                                                      | Device: 119 g (toothbrush + 2 brush heads) |
| Battery:<br>Capacity<br>Nominal voltage<br>Type designation | 600 mAh<br>3.7 V<br>Lithium-ion battery    |
| Power supply                                                | 5.0 V $\overline{\text{---}}$ 1.0 A, 5.0 W |

Subject to technical changes

11. GUARANTEE

Further information on the guarantee and guarantee conditions can be found in the guarantee leaflet supplied.

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Beurer GmbH • Söflinger Str. 218
BEURER UK Ltd., Suite 16,